

EXP PRE-COLUMN FILTER

The EXP Pre-Column Hand-Tight Filter is ideal for protecting small particle size UHPLC columns. Such techniques analyze samples in the most demanding applications, which can decrease the life of these expensive columns. EXP Pre-Column Filters help extend that life and protect your column investment without sacrificing performance. Free-turn® architecture allows the user to **change cartridges by hand** without breaking fluid connections on the holder inlet/outlet.

EXP Pre-Column Filter holders are available with Titanium Hybrid ferrules for easy direct connection to any 10-32 port. The filter holder comes as a complete package including fittings to provide repeated zero-dead-volume column connections.


ADVANTAGES OF EXP PRE-COLUMN FILTERS

- Approved for use at 20,000 psi
- Tested to 30,000 psi
- Depth filtration maximizes capacity
- Auto-adjusting ZDV column connection
- Custom volumes/porosities available
- Hand-tight filter replacement - NO TOOLS!
- Low-volume, low-dispersion cartridges
- Hardened stainless steel end cap eliminates galling

simple.
innovative.
easy to use.

UHPLC

www.optimizech.com/exp

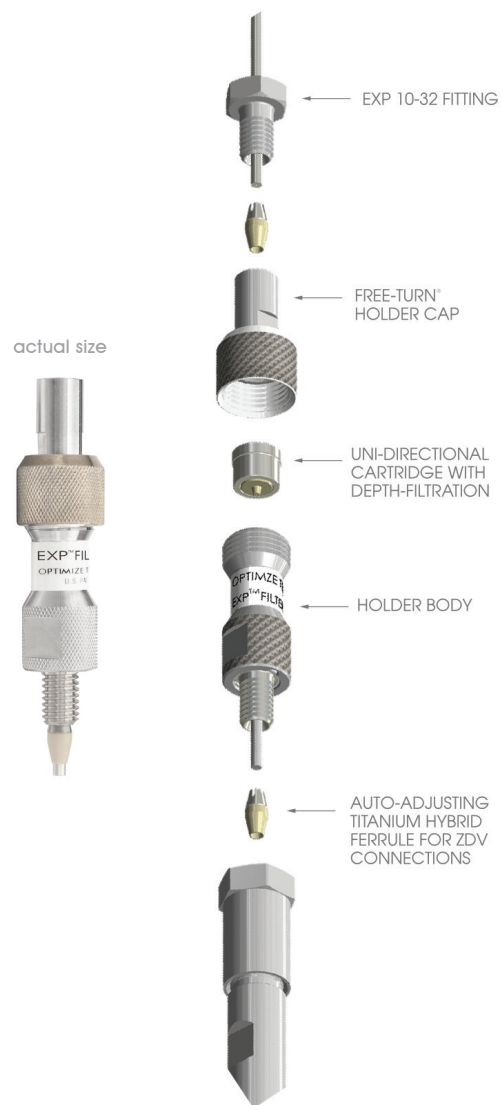

EXP PRE-COLUMN FILTER HOLDER

15-04-03837

EXP Pre-Column Filter Holder
(fittings included)



EXP PRE-COLUMN REPLACEMENT FILTER CARTRIDGES

15-04-03093	0.5 µm, EXP Filter Cartridge, 5/Pk
15-04-03096	0.2 µm, EXP Filter Cartridge, 5/Pk
15-04-03094	0.5 µm, EXP Filter Cartridge, 10/Pk
15-04-03097	0.2 µm, EXP Filter Cartridge, 10/Pk